

HY-BOND GLASIONOMER ***CX-Smart***

Smart Self-Adhesive Cementation



HY-BOND GLASIONOMER **CX-Smart**

Quick **Easy**
Predictable

Hy-bond Glaslonomer CX-Smart is a truly biocompatible, quick-set, contemporary luting glass ionomer that guarantees reliable and lasting restorations.

Based on Shofu's time-tested, self-adhesive luting technology, **Hy-bond CX-Smart** comprises of Shofu's proprietary additive - **Hy-agent**, "a tannic acid-strontium fluoride-zinc fluoride complex" that is known to **promote re-mineralization**, **increase acid resistance** and **reduce solubility** with anti-bacterial properties while virtually **eliminating post-op hypersensitivity**.

BENEFITS FOR AN EFFICIENT CLINICAL WORKFLOW



- *Exceptional biocompatibility with favourable tissue response*
- *Chemically bonds to the tooth structure for excellent marginal seal*
- *Quick-set with remarkable moisture tolerance*
- *No etch or bond for a quick, easy and simple application – ideal for paediatric, geriatric and special needs patients.*
- *Resin-free, suitable for patients with resin allergy*
- *Optimal thixotropicity and even film thickness*
- *Durable physical properties, high acid resistance and low oral solubility*
- *Hy-additive ensures virtually no post-operative sensitivity & favourable pulp response*
- *Cost-effective and efficient*

IDEAL FOR CEMENTATION OF

	Metal	Porcelain Fused to Metal	Zirconia
Inlays / Onlays	✓		
Crown & Bridges	✓	✓	✓
Post & Core	✓		
Orthodontic Appliances	✓		

DURABLE RESTORATIONS

Hy-bond CX-Smart exhibits a robust matrix that resists crack propagation and greater acid resistance for long-term stability.

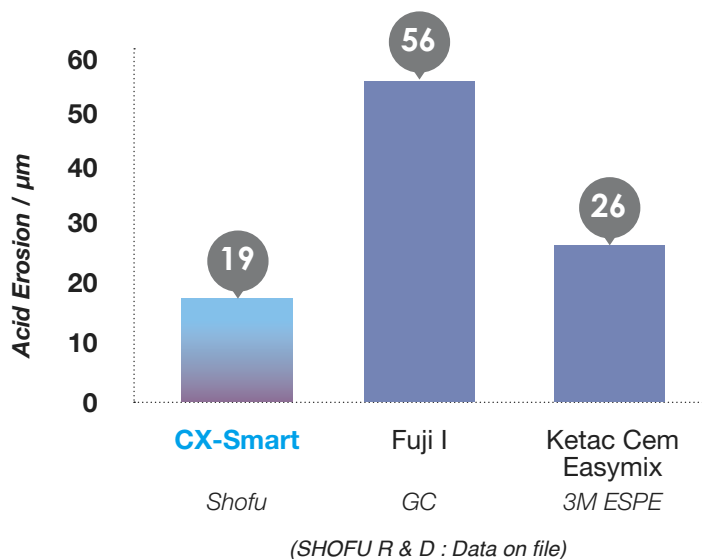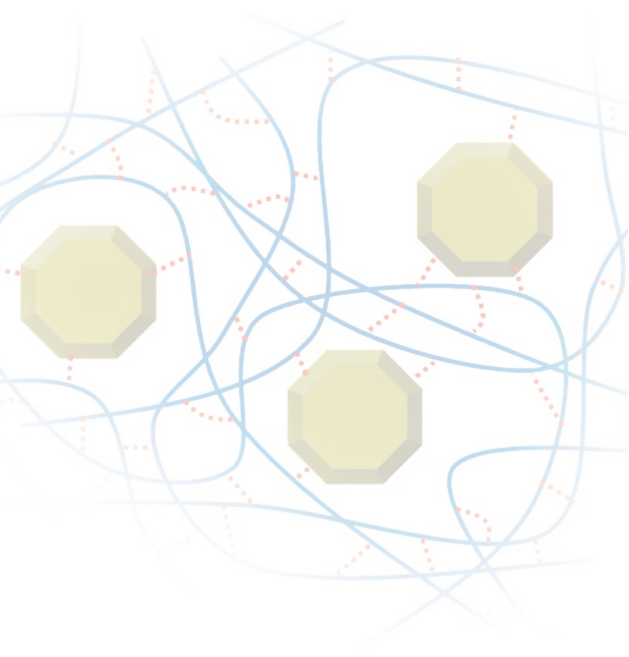

SMARTER & EFFICIENT WORK FLOW

A polyvalent acrylic-tricarboxylic acid co-polymer formulation imparts ideal thixotropicity, so the cement easily flows under pressure to a uniform film thickness for a perfect fit.



Photographs courtesy of Dr. Rintaro Terata



Hy-bond CX-Smart exhibits a sharp clinical set reducing chair time and effortless removal of excess cement for a clean finish and greater patient comfort.

Clinical Time

Mixing time (23 °C / 73 °F)	30"	—
Working time (23 °C / 73 °F)	2' 00"	From start of mixing
Setting time (37 °C / 99 °F)	2' 30"	From end of mixing
Removal of set cement (rubbery consistency)	1' 30"	After seating restoration

*Note: Setting time after end of mixing stipulated in ISO 9917-1 is 2'30"


Order Information



Hy-Bond Glaslonomer CX-Smart 1-1 Set
Powder 45gm / Liquid 30gm (25ml)
[PN C0017]

Hy-Bond Glaslonomer CX-Smart Mini Set
Powder 15gm / Liquid 10gm (8ml)
[PN C0018]

Related Products

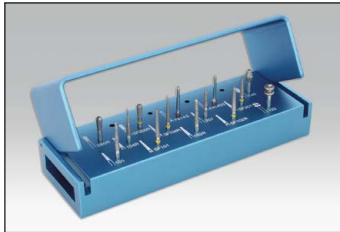

Tooth Preparation :
Crown & Bridge Preparation Kit
[PN F270]



Chair-side Adjustment :
Porcelain Adjustment Kit HP
[PN 0301]



Cement Removal :
Silicone One Gloss Midi
[PN 0171]



Intra-oral Adjustment :
Enamel Adjustment Kit CA
[PN 0307]

Find us at



SHOFU INC. 11 Kamitakamatsu-cho, Fukuine, Higashiyama-ku, Kyoto 605 - 0983, Japan (www.shofu.co.jp)
 SHOFU DENTAL ASIA-PACIFIC PTE. LTD. 10 Science Park Road, #03-12 The Alpha, Science Park II, Singapore 117684 (www.shofu.com.sg)
 SHOFU DENTAL TRADING (SHANGHAI) CO. LTD. No.645 Jiye Road, Sheshan Industry Park Songjiang Shanghai 201602, China (www.shofu.com.cn)
 SHOFU DENTAL CORPORATION 1225 Stone Drive, San Marcos, CA 92078-4059, USA (www.shofu.com)
 SHOFU DENTAL GmbH Am Brüll 17, 40878 Ratingen, Germany (www.shofu.de)
 SHOFU UK. Riverside House, River Lawn Road, Tonbridge, Kent TN9 1EP, United Kingdom (www.shofu.co.uk)