




Sturdy Articulator Designed for Guaranteed Performance

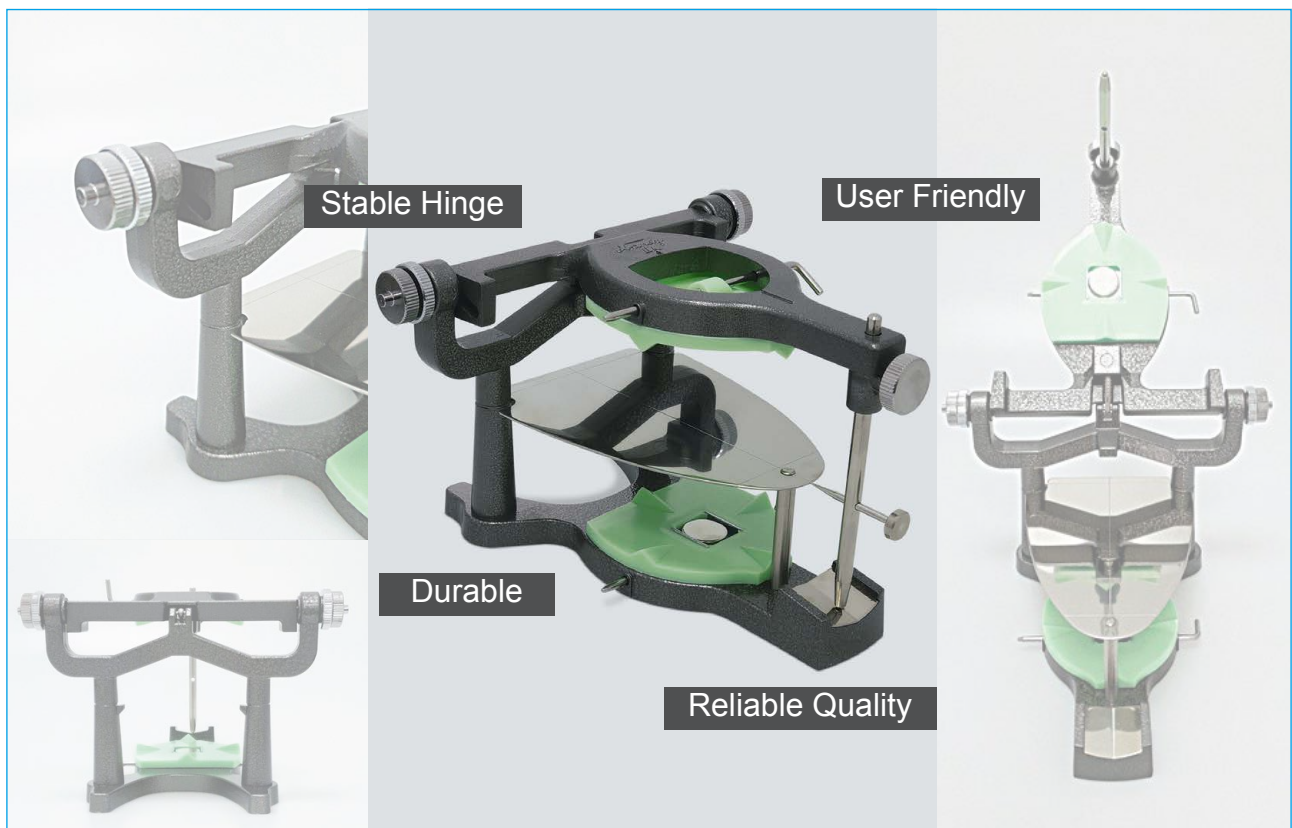

Stable Hinge • Reliable Quality • Durable • User Friendly

Trusted for Decades to Create Prosthetic work with Precision



## ***Make Handy IIA Articulator your Preferred Choice Today !***


### **Features**

- Easy to use, reliable, and durable
- Able to work with more than one model at a time with 2 quick removable magnetic split cast mounting plates
- Permits hinge movement when the centric screw is locked and permits condylar movement when the centric screw is loosened
- Stable Hinge provides accurate condylar movement each time
- Upper jaw can be inclined to 45° for convenient mounting of the model
- The horizontal position can be maintained even when the Incisor Rod is removed

### **Specifications**

Incline of the condylar path	25°
Bonwill triangle	25°
Distance between condyles	10.5cm
Distance between condyles & central incisor of lower jaw	11.2cm
Distance between the upper and lower jaw	8cm

***Handy IIA Articulator has been used and recommended by Dental Schools and Dental Laboratories around the world ...***

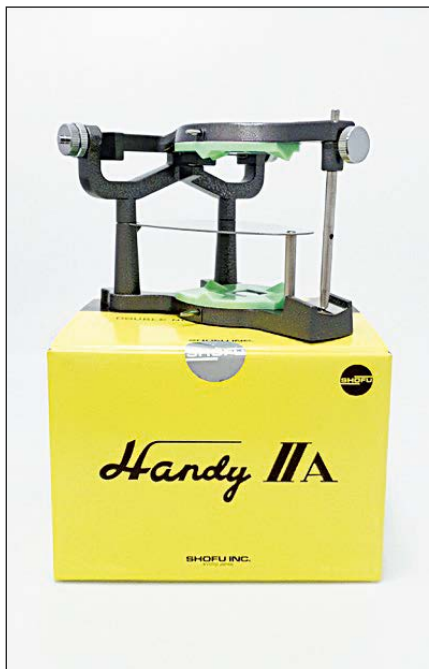

- Individually packed [PN 5041]



**SHOFU DENTAL ASIA-PACIFIC PTE. LTD.**  
[www.shofu.com.sg](http://www.shofu.com.sg)