

CLICK / SCAN




FOR MORE INFO

LAB CARBIDE

DURABLE & EFFICIENT

A reliable range of carbide burs in different shapes and flute designs for fast chair-side or bench side adjustment of various dental substrates.

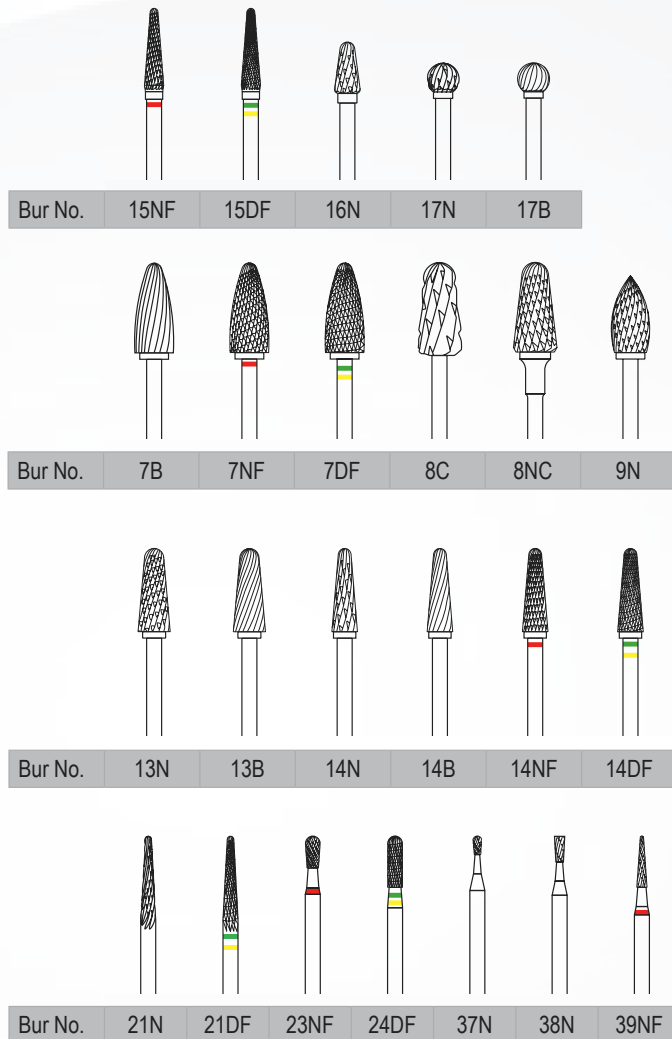

Available in 31 shapes
Individually packed



APPLICATIONS

Adjustment and finishing of

- Acrylic Resins
- Nickel-Chrome Cobalt Alloys
- Chrome-Nickel Alloys
- Precious Alloys

FEATURES

- ✓ **Enhanced Cutting Efficiency**
High-grade tungsten carbide burs with remarkable durability
- ✓ **Uniform Concentricity**
Minimal vibration and precise trimming / adjustment of either acrylic resin or cast metal
- ✓ **6 different types**
 - Coarse (C)
 - Normal Coarse (NC)
 - Normal Cross Cut (N)
 - Basic Plain Cut (B)
 - Normal Fine Cut (NF)
 - Diamond Fine Cut (DF)

DIRECTION FOR USE

1. C – Gross reduction of plaster
2. NC , N , B – For trimming and initial finishing
3. NF & DF – For anatomical contouring and finishing

RELATED PRODUCTS

- AcryPoints
- Acrylic Contouring & Finishing Kit (REF 0329)

